

Quick reference guide for SHKE8000-8CE

Step 1: Pull out the platform (fig.1)

- Open the chamber door.
- Grasp the center handle, pull the inner release bar upward, and then press the handle downward until it clicks into a horizontal position.
- Pull outward on the platform.

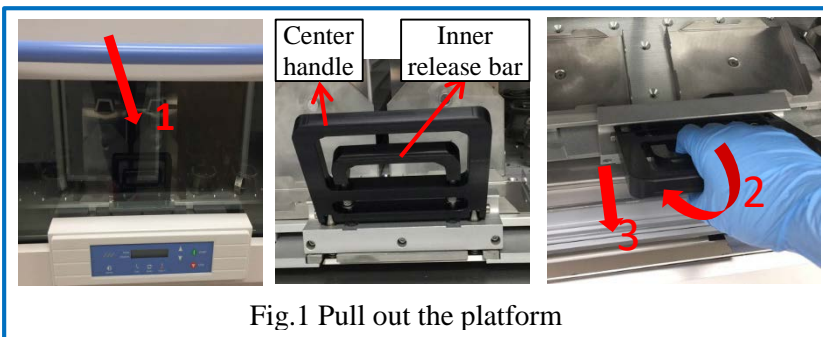


Fig.1 Pull out the platform

Step2: Load your samples, making sure all flasks and test tube racks are firmly seated in the clips.

Caution: *Don't operate this shaker with flammable materials.*

Step3: Push the platform all the way into the chamber. Pull up on the inner release bar and press the center handle upward into a vertical position. Close the chamber door.

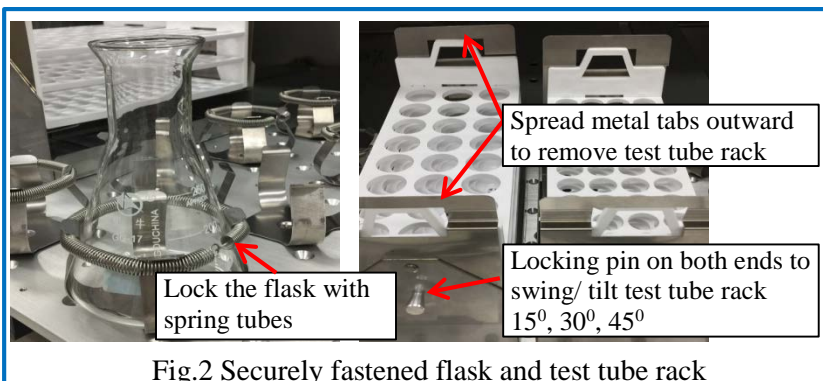


Fig.2 Securely fastened flask and test tube rack

Step4: Power on the shaker.

Step 5: Control panel operation

- Press / / , then press or to set time/speed/ run temperature;
- Press twice, then press or to set the hold temperature;



Fig.3 Control panel

Step 6: Press to start.

Caution: *It may take up to one minute to bring the platform up to full speed. Never leave the shaker unattended while starting.*

Step 7: Press to stop when finished if the time mode is hold.